

DESCRIPTION

Target:	OLR1
Target aliases:	CLEC8A, SCARE1, LOX-1, HLOX-1, LOX1, LOXIN, SLOX1
Fc isotype:	IgG2a
Membrane proteome specificity:	Monospecific for 6,000 membrane proteins tested
Species reactivity:	Human (Others untested)
Extracellular epitope:	Yes
Fc modifications:	C-terminal Avitag ¹ , disabled Fc-γ receptor binding ²
Source:	Recombinant CHO expression; purified by Protein A chromatography
Formulation:	Endotoxin Free PBS pH 7.4, sterile-filtered
Concentration:	1 mg/ml

1. A peptide tag that can be biotinylated in vitro using the biotin ligase enzyme (BirA).
2. Mutated Fc-γ receptor binding site to minimize non-specific antibody binding.

OLR1 TARGET INFORMATION


OLR1 is a single-pass transmembrane protein and low density lipoprotein receptor with a secreted form. OLR1 is part of the C-type lectin superfamily. OLR1 is responsible for binding, internalizing, and degrading oxidatively modified low density lipoprotein (oxLDL) and possibly regulates Fas-induced apoptosis. Mutations in OLR1 are associated with atherosclerosis and myocardial infarction risk and may be related to Alzheimer's disease risk. (NCBI Gene: 4973, UniProtKB/Swiss-Prot: P78380) Other names: CLEC8A, SCARE1, LOX-1, HLOX-1, LOX1, LOXIN, SLOX1

SHIPPING AND STORAGE

Shipping:	Shipped at ambient temperature. Store at 4°C.
Stability & Storage:	Stable for 12 months from date of receipt when stored at 4°C. Avoid repeated freeze-thaw cycles.

VALIDATION DATA


Membrane Proteome Specificity

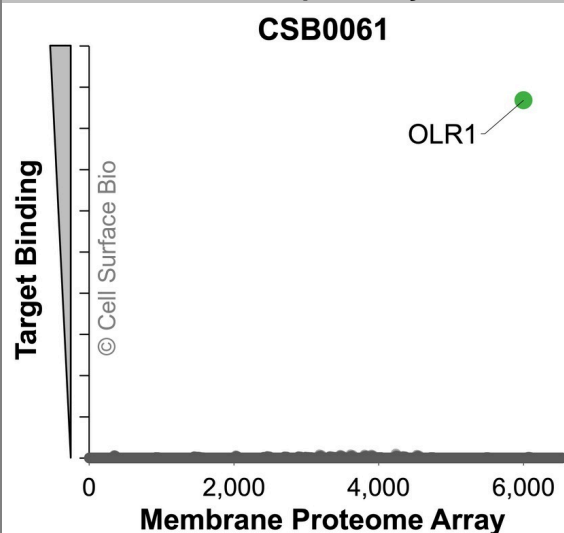

The specificity of OLR1 Monoclonal Antibody (CSB0061) was tested on the Membrane Proteome Array™ and shown to be specific for OLR1.


The Membrane Proteome Array™ contains 6,000 different human membrane proteins, each expressed in unfixed human cells to ensure native conformation and post-translational modifications. The Membrane Proteome Array™ represents the industry standard for determining the binding specificity of antibodies and other protein ligands.

Applications	Conditions	Recommended concentration
Flow Cytometry, Extracellular	Live, Unpermeabilized	1 µg/ml

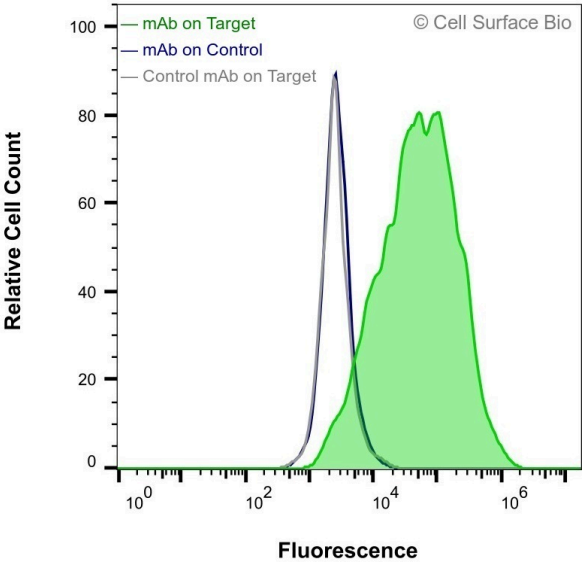


HEK-293F cells transiently transfected with OLR1 were stained with OLR1 Monoclonal Antibody (CSB0061) (green) or isotype control antibody (gray), followed by AlexaFluor 647- conjugated anti- IgG secondary antibody. HEK-293F cells transiently transfected with an empty control vector were also stained with OLR1 Monoclonal Antibody (CSB0061) (blue).

Applications	Conditions	Recommended concentration
Immunofluorescence, Extracellular	Fixed 4% paraformaldehyde	1 µg/ml

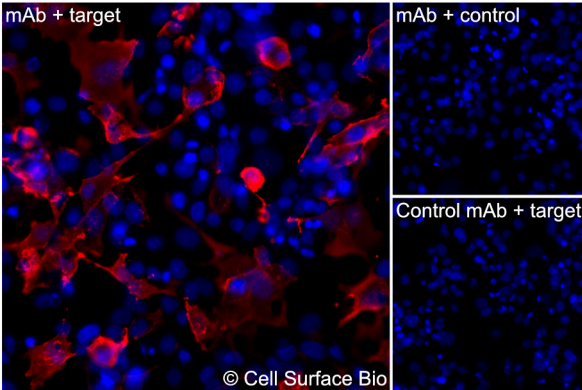

(Left) COS-7 cells transiently transfected OLR1 were stained with OLR1 Monoclonal Antibody (CSB0061) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with OLR1 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with OLR1 and stained with control MAb.

Applications	Conditions	Recommended concentration
Immunofluorescence, Intracellular	Fixed 4% paraformaldehyde, Permeabilized 0.1% Triton X-100	1 µg/ml

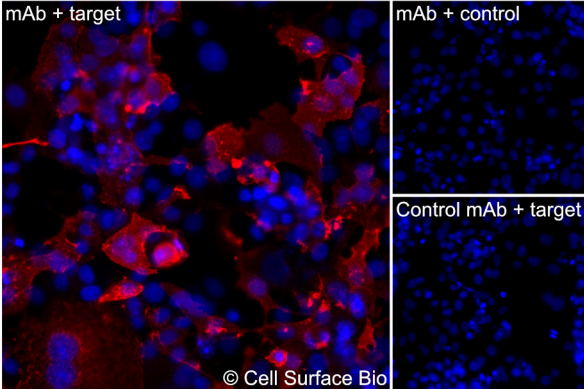

(Left) COS-7 cells transiently transfected OLR1 were permeabilized and stained with OLR1 Monoclonal Antibody (CSB0061) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with OLR1 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with {{target}} and stained with control MAb.