

DESCRIPTION

Target:	Chikungunya Virus E1/E2
Target aliases:	
Fc isotype:	IgG2a
Membrane proteome specificity:	Monospecific for 6,000 membrane proteins tested
Species reactivity:	Virus
Extracellular epitope:	Yes
Fc modifications:	C-terminal Avitag ¹ , disabled Fc- γ receptor binding ²
Source:	Recombinant CHO expression; purified by Protein A chromatography
Formulation:	Endotoxin Free PBS pH 7.4, sterile-filtered
Concentration:	1 mg/ml


1. A peptide tag that can be biotinylated in vitro using the biotin ligase enzyme (BirA).
2. Mutated Fc- γ receptor binding site to minimize non-specific antibody binding.

Chikungunya Virus E1/E2 TARGET INFORMATION

Chikungunya virus (CHIKV) is spread through bites from infected mosquitoes and can cause fever, joint pain and swelling, headache, muscle pain, or a rash (<https://www.cdc.gov/chikungunya/index.html>).

SHIPPING AND STORAGE


Shipping:	Shipped at ambient temperature. Store at 4°C.
Stability & Storage:	Stable for 12 months from date of receipt when stored at 4°C. Avoid repeated freeze-thaw cycles.

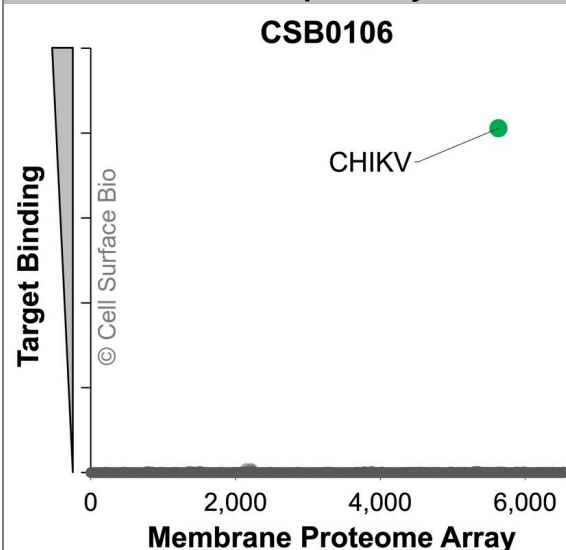VALIDATION DATA
Membrane Proteome Specificity


The specificity of Chikungunya Virus E1/E2 Monoclonal Antibody (CSB0106) was tested on the Membrane Proteome Array™ and shown to be specific for Chikungunya Virus E1/E2.


The Membrane Proteome Array™ contains 6,000 different human membrane proteins, each expressed in unfixed human cells to ensure native conformation and post-translational modifications. The Membrane Proteome Array™ represents the industry standard for determining the binding specificity of antibodies and other protein ligands.

Applications	Conditions	Recommended concentration
Flow Cytometry, Extracellular	Live, Unpermeabilized	1 µg/ml

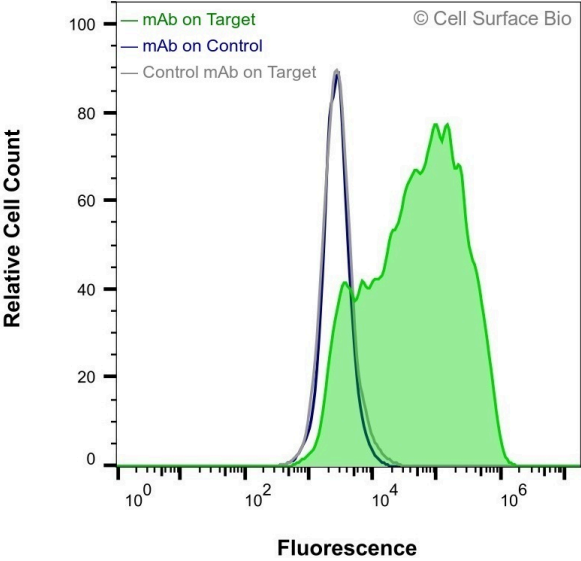


HEK-293F cells transiently transfected with Chikungunya Virus E1/E2 were stained with Chikungunya Virus E1/E2 Monoclonal Antibody (CSB0106)(green) or isotype control antibody (gray), followed by AlexaFluor 647- conjugated anti- IgG secondary antibody. HEK-293F cells transiently transfected with an empty control vector were also stained with Chikungunya Virus E1/E2 Monoclonal Antibody (CSB0106)(blue).

Applications	Conditions	Recommended concentration
Immunofluorescence, Extracellular	Fixed 4% paraformaldehyde	1 µg/ml

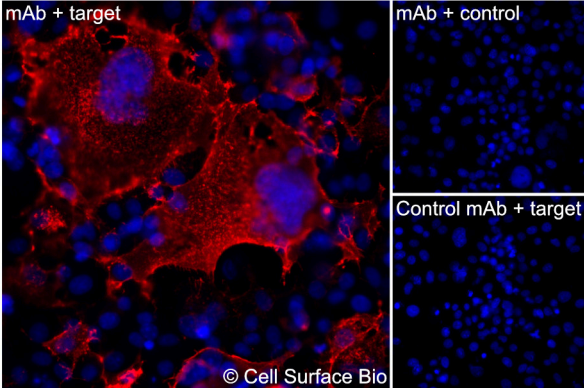

(Left) COS-7 cells transiently transfected Chikungunya Virus E1/E2 were stained with Chikungunya Virus E1/E2 Monoclonal Antibody (CSB0106) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with Chikungunya Virus E1/E2 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with Chikungunya Virus E1/E2 and stained with control MAb.

Applications	Conditions	Recommended concentration
Immunofluorescence, Intracellular	Fixed 4% paraformaldehyde, Permeabilized 0.1% Triton X-100	1 µg/ml

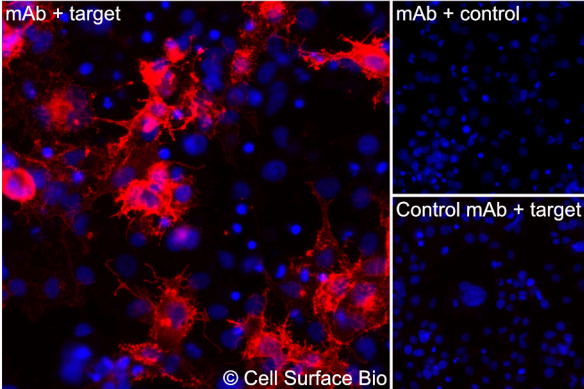

(Left) COS-7 cells transiently transfected Chikungunya Virus E1/E2 were permeabilized and stained with Chikungunya Virus E1/E2 Monoclonal Antibody (CSB0106) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with Chikungunya Virus E1/E2 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with {{target}} and stained with control MAb.