

DESCRIPTION

Target:	CD20
Target aliases:	MS4A1, Bp35, FMC7, B1, LEU-16, CVID5, S7
Fc isotype:	IgG2a
Membrane proteome specificity:	Monospecific for 6,000 membrane proteins tested
Species reactivity:	Human (Others untested)
Extracellular epitope:	Yes
Fc modifications:	C-terminal Avitag ¹ , disabled Fc- γ receptor binding ²
Source:	Recombinant CHO expression; purified by Protein A chromatography
Formulation:	Endotoxin Free PBS pH 7.4, sterile-filtered
Concentration:	1 mg/ml


1. A peptide tag that can be biotinylated in vitro using the biotin ligase enzyme (BirA).
2. Mutated Fc- γ receptor binding site to minimize non-specific antibody binding.

CD20 TARGET INFORMATION


CD20 is a B-lymphocyte membrane protein that is part of the membrane-spanning 4A gene family. CD20 is involved in regulating calcium influx, which is important for the development and differentiation of B-cells into plasma cells. (NCBI Gene: 931, UniProtKB/Swiss-Prot: P11836) Other names: MS4A1, Bp35, FMC7, B1, LEU-16, CVID5, S7

SHIPPING AND STORAGE

Shipping:	Shipped at ambient temperature. Store at 4°C.
Stability & Storage:	Stable for 12 months from date of receipt when stored at 4°C. Avoid repeated freeze-thaw cycles.

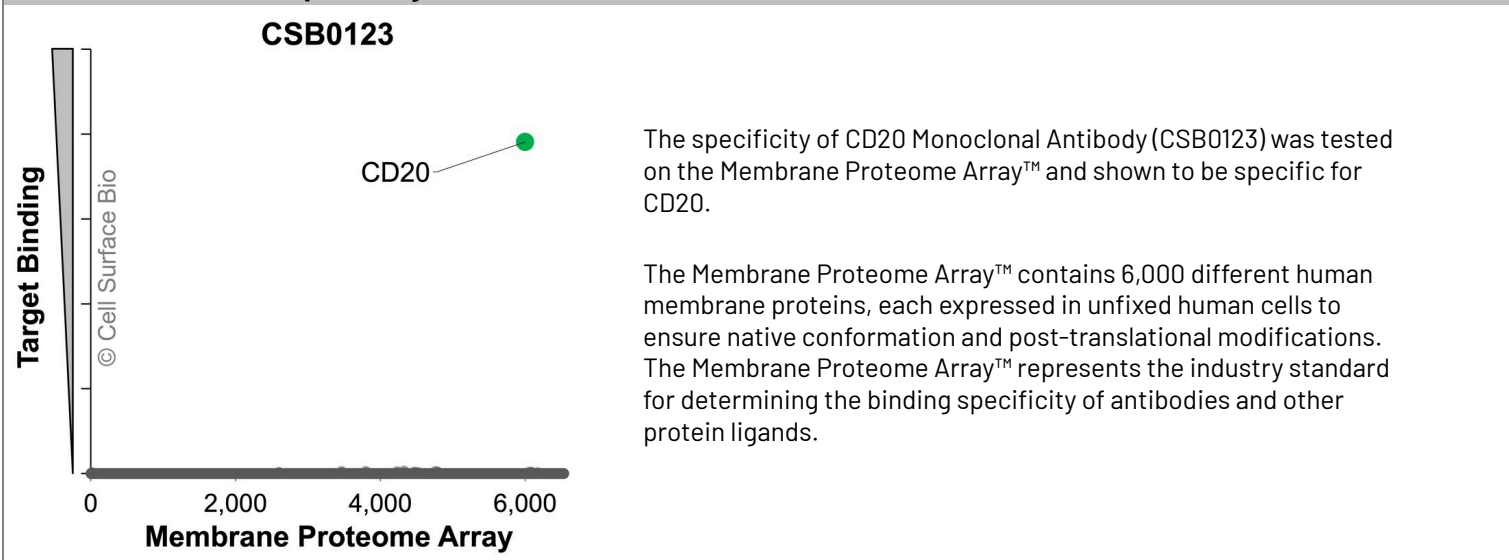
VALIDATION DATA
Membrane Proteome Specificity


Applications	Conditions	Recommended concentration
Flow Cytometry, Extracellular	Live, Unpermeabilized	1 µg/ml

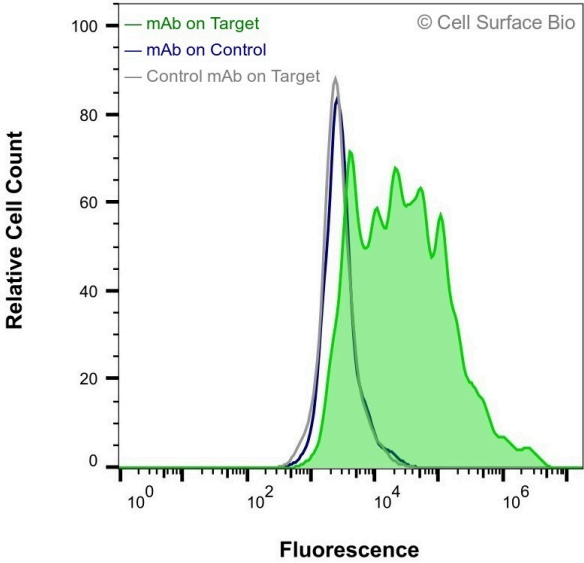

HEK-293F cells transiently transfected with CD20 were stained with CD20 Monoclonal Antibody (CSB0123) (green) or isotype control antibody (gray), followed by AlexaFluor 647- conjugated anti- IgG secondary antibody. HEK-293F cells transiently transfected with an empty control vector were also stained with CD20 Monoclonal Antibody (CSB0123) (blue).

Applications	Conditions	Recommended concentration
Immunofluorescence, Extracellular	Fixed 4% paraformaldehyde	1 µg/ml

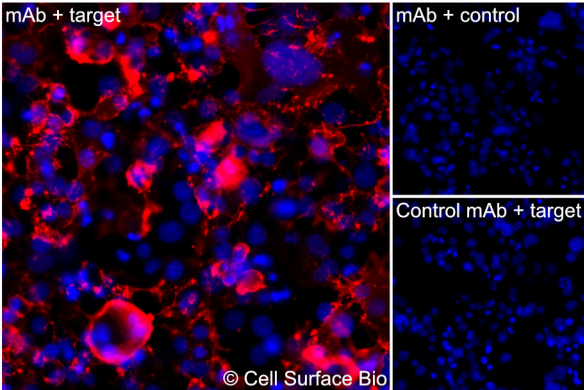

(Left) COS-7 cells transiently transfected CD20 were stained with CD20 Monoclonal Antibody (CSB0123) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with CD20 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with CD20 and stained with control MAb.

Applications


Immunofluorescence, Intracellular

Conditions

Fixed 4% paraformaldehyde,
Permeabilized 0.1% Triton X-100

Recommended concentration

1 µg/ml

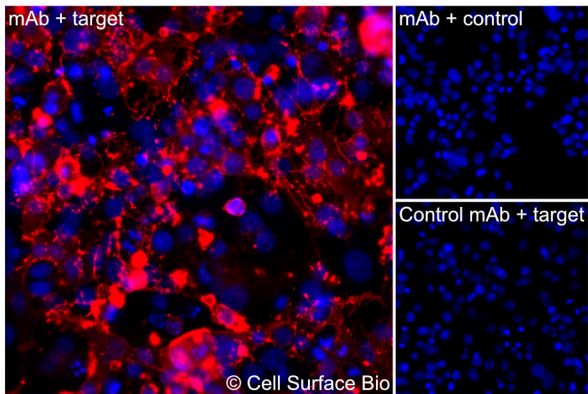

(Left) COS-7 cells transiently transfected CD20 were permeabilized and stained with CD20 Monoclonal Antibody (CSB0123) followed by AlexaFluor 647 anti- IgG secondary antibody (red) and DAPI (blue). (Top right) COS-7 cells transiently transfected with an empty control vector stained with CD20 Monoclonal Antibody. (Bottom right) Isotype control: COS-7 cells transfected with {{target}} and stained with control MAb.